



incubators  
series 636 PLUS  
NAHITA

The NAHITA incubators Serie 636 plus are used in a lot of process in laboratories of cosmetics, microbiology, food analysis, biotechnology, botany and biochemistry, also is used in some applications in serum and plasma test, conservation or incubation of samples and the maintenance of cell cultures and microorganisms.

Made of painted steel and kiln-dried, have a stainless steel interior for an easy clean and a door with double glazed window that lets you control the evolution of the samples at any times.

They have digital regulation of temperature and time so that allow a precise control of these processes.


## features and design of the equipment

Modern and functional design destined to extract the maximum output to your equipment:

- Exterior metal furniture cover with epoxy and interior of stainless steel, highly resistant to acids and corrosives, very easy to clean.
- Door with double glazed window for external control, no temperature variations.
- Supplied with trays (2 units, height adjustable), which provides great flexibility of use allowing to introduce bottles and containers of different sizes.



# incubators series 636 PLUS: equipment parts

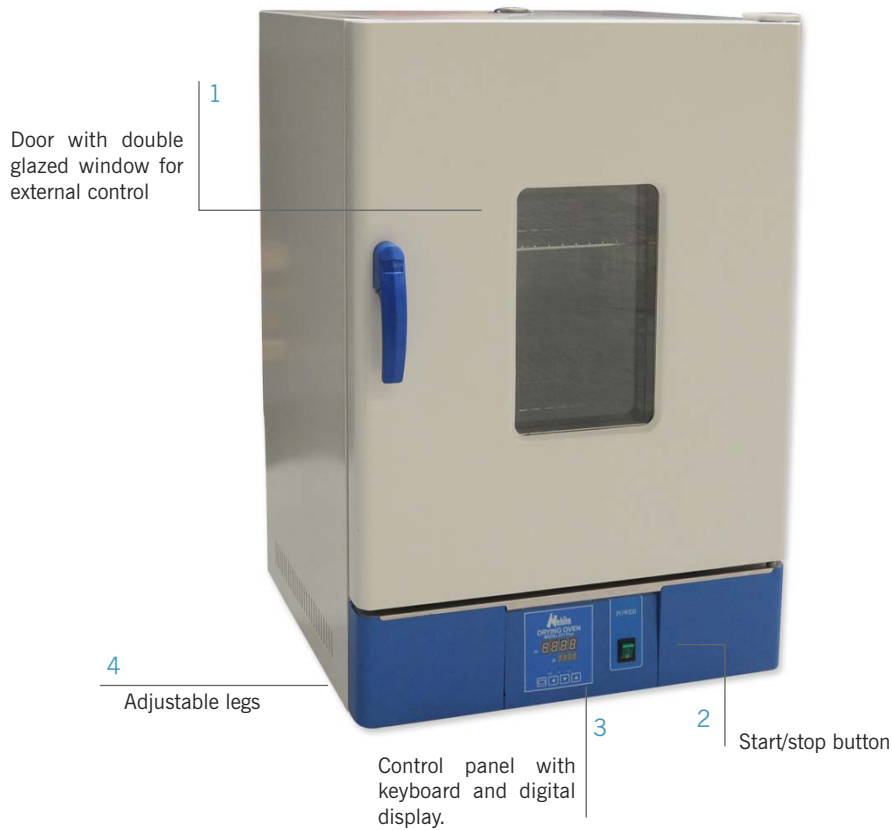

incubatorsNAHITA

# incubators series 636 PLUS: back of the equipment



# incubators series 636 PLUS: front panel

- [01] Temperature display
- [02] Display user-selected temperature
- [03] Alarm light LED
- [04] Light AutoFit LED
- [05] Operating light LED
- [06] Set-up button LED
- [07] Cursor button
- [08][09] Time and Temperature button
- [10] Switch ON / OFF



## technical specifications

### Specifications

- Air circulation system by natural convection
- Control panel with keyboard and digital display.
- Adjustable legs for easy installation.
- Protection system overheat with visual and audible alarm.

## incubators series 636 PLUS: capacidades y modelos

[01] Code 50636031  
Capacity: 30 L



[02] Code 50636041  
Capacity: 45 L



[03] Code 50636061  
Capacity: 65 L



[04] Code 50636131  
Capacity: 125 L



## table technical specifications

Code	50636031	50636041	50636061	50636131
Model	636/3 Plus	636/45 Plus	636/75 Plus	636/125 Plus
Capacity	30 L	45 L	65 L	125 L
Max. Temperature	RT+5-80 °C			
Resolution	0.1 °C			
Temperature fluctuation	± 0.5 °C			
Power	300 W	300 W	400 W	400 W
Dimensions (cm)				
Exterior	46x51x69	50x55x73	55x55x84	64x68x91
Interior	31x31x31	35x35x35	40x36x45	50x45x55
Weight	39 Kg	42 Kg	53 Kg	58 Kg
Fuses	15 A			
Power supply	220 V, 50 Hz ± 10%			



**AUXILAB S.L.**

Material de laboratorio  
Laboratory supplies

Polígono Morea Norte  
Calle D Nº 6 Beriain [Navarra]  
C.P. 31191  
Tel: +34 948 310 513  
[www.auxilab.es](http://www.auxilab.es)